

Exclusively Distributed by



PRISM MEDICAL UK


Gantry Systems

FSG200 Free Standing Gantry

The FSG200 Free Standing Gantry system is designed to provide a safe, stable hoist track in places where a fixed ceiling track is not available or portable solution is required.

When used in conjunction with a portable track hoist, such as the CP200, the system is perfect for home care, allowing users to extend the time they remain in their own home by years. It's a cost effective solution that is easy to transport and simple erect in minutes by either 1 or 2 people without the need for tools.

With complete safety, dignity and comfort for the client and carer the FSG200 Free Standing Gantry requires no structural alterations to the room or ceiling being effectively used to lift and transfer individuals from a bed, chair or similar fixture.

Key Features

1. Safe Working Load (SWL) - 200kg
2. Compatible with CP200 Portable Track Hoist
3. High-strength lightweight aluminium
4. Adjustable in length from 1.9 - 2.85m, and in height from 2.1 - 2.5m
5. Stable 0.87m 'feet'
6. Built-in easy-slide trolley
7. Floor level adjusters
8. Star Handle securing screws
9. Safety latch to prevent uprights retracting
10. Securing Bolts and lynch pins



1 Built in Easy-Slide Hoist Trolley Set




2 Floor levelling adjusters on each foot



3 Safety latches to prevent uprights from accidentally retracting



4 Star-Handle securing screws on track section and feet

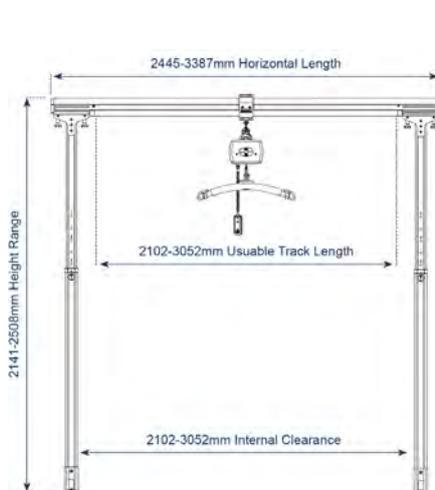

Code	Description	Useable Track Length	Width	Height	Foot Length	Unit Weight	Height Adjustments	SWL
LIC352105	FSG200 Free Standing Gantry	1900 - 2825mm	2445 - 3387mm	2141 - 2508mm	870mm	47.5kg	75mm	200kg

FSG200 Free Standing Gantry is designed for use in conjunction with:



P200 Ceiling Hoist
(LIC350110 / LIC350120)



CP200 Ceiling Hoist
(LIC350105)



✓ FSG200 Gantry is registered in the ARTG, AUST R ARTG 286310

Pressure Fit System

Pressure Fit System

The unique and portable Pressure Fit System is light and compact – manufactured from high strength aluminium, each component part weighs no more than 12kg. The system is easy to transport and assemble, even where space is limited and other systems cannot be used.

Minimal installation is required because the system works by exerting equal pressure between the floor and ceiling to stabilise the unit.

It can be used to hoist between beds, chairs and toilets/ commodes, either as a temporary or permanent solution.

Key Features


1. Safe Working Load (SWL) – 200kg
2. Works with ceilings from 2133mm to 2743mm
3. Span adjustable from 2100mm to 3000mm
4. Adjustable in length from 1.85 - 2.75m, and in height from 2 - 2.75m
5. Easy to read pressure gauge indicator ensures safe installation
6. Can be used with all Freeway portable ceiling lifts



1 Built in Easy-Slide Hoist Trolley Set with QRS connector loop




2 Ratchet system pressure fits uprights between floor and ceiling

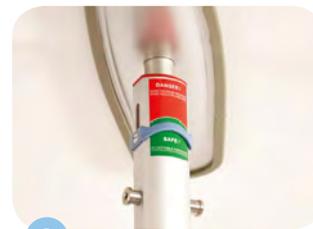

3 Pressure fit safety gauge ensures adequate pressure has been applied



4 Bubble level to ensure uprights are vertical

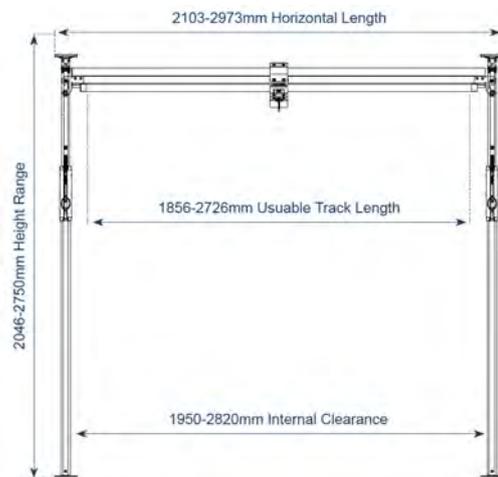

Code	Description	Useable Track Length	Width	Height	Unit Weight	SWL
LIC352120	Pressure Fit System	1856 - 2726mm	2103 - 2973mm	2046 - 2750mm	27.1kg	200kg

Pressure Fit System is designed for use in conjunction with:



P200 Ceiling Hoist
(LIC350110 / LIC350120)



CP200 Ceiling Hoist
(LIC350105)



Pressure Fit System is registered in the ARTG, AUSTRALIA ARTG 286310

Exclusively Distributed by

 aidacare

1300 133 120

www.aidacare.com.au



PRISM MEDICAL UK